

# RETRACTABED®

MED MIZER®

The Ultimate Choice  
For Space Saving Design

Med-Mizer SS Series  
**RETRACTABED**  
Resident Safety, Convenience  
and Independence



### Fall Risk:

Non-Retractable Bed

### Fall Prevention:

RetractaBed

## Resident Safety

In a non-retractable bed, the residents place themselves in a dangerous position of stretching out over the bed to answer the telephone, reach for glasses, or for a glass of water. Thus, normal resident movements are a constant threat to resident safety. The RetractableBed® will significantly reduce these dangers to residents.

## Resident Convenience

When the one button Auto Contour position is actuated, the head section of the bed hugs the wall, thus keeping the resident in proximity to the bedside cabinet and overhead lighting.

## Resident Independence

The independence offered to the resident with the RetractableBed® gives the resident complete control of the resident-room environment. This sophisticated retraction feature is the solution to present and future facility personnel shortages.



## Exclusive Ergonomically Correct Med-Support Arm

**Arm Positions:** The low position allows unobstructed access to residents or for stretcher transfer.

The normal guarded position provides protection without restraints and a pleasing and unprecedented arm chair environment when contoured.

The third position with the arm high and towards the head end provides the resident with an extremely easy and stable means of entering and exiting the bed.

**Easy Access Support Arm Controls:** Bed Positioning controls are conveniently located on both support arms. The Clinical-Contour, or head & knee functions, are located on the inside of both support arms for easy access by the patient or resident. The Hi/Low Control function is located on the outside of both support arms for easy caregiver access.

**Point Of Care:** Fixed/Embedded controls and a three position med-support arm allow the staff and resident to be in a position to maximize the safety of their environment. This exclusive feature to the RetractaBed also gives residents the capability to move freely in and out of bed thus creating an assistance device not a restraint.

*The Med-Support Arm complies with all Entrapment zones*

## Retractabed Clinical - Contour™ Feature

By activating a single control button, the head section and knee sections simultaneously create a comfortable 12 to 15 degree knee break. This auto-contour feature prevents the resident from gravitating toward the foot of the bed when the head section is elevated or lowered. This eliminates many problems nursing homes face today: skin shear leading to skin ulcers, repositioning residents in bed leading to skin shearing and possible back injuries to staff. This feature also relieves undue pressure behind the knees, thus assuring proper circulation to the lower extremities.


**No Clinical-Contour™:** As a resident's position changes, 5 to 7 inches of gravitation to the footboard can occur, a leading cause of skin shearing in the lower back and ankles. Stress is also placed on the patella, hindering blood flow to the lower extremities.



**Clinical-Contour™:** Position of resident stays in relatively the same position in the bed, reducing skin shearing which may cause skin ulcers. This helps relieve pressure behind the knees, providing better blood flow to the lower extremities. Staff risk of injury due to repositioning residents is lowered with Clinical-Contour™.



## Accessories

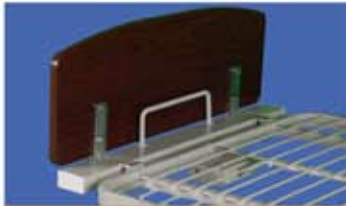


4" Bed Extender



5" Caster Upgrade




Fixed Assist Bar

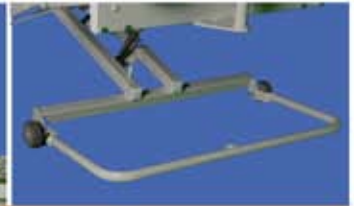

Floor Bumper Guard



Med Support Arm Pad

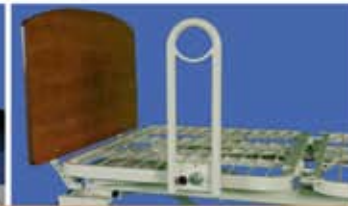

Pivoting Assist Bars




Trapeze Assembly




Transport Dolly

## Features

Pendant Control  
(Optional, not shown)

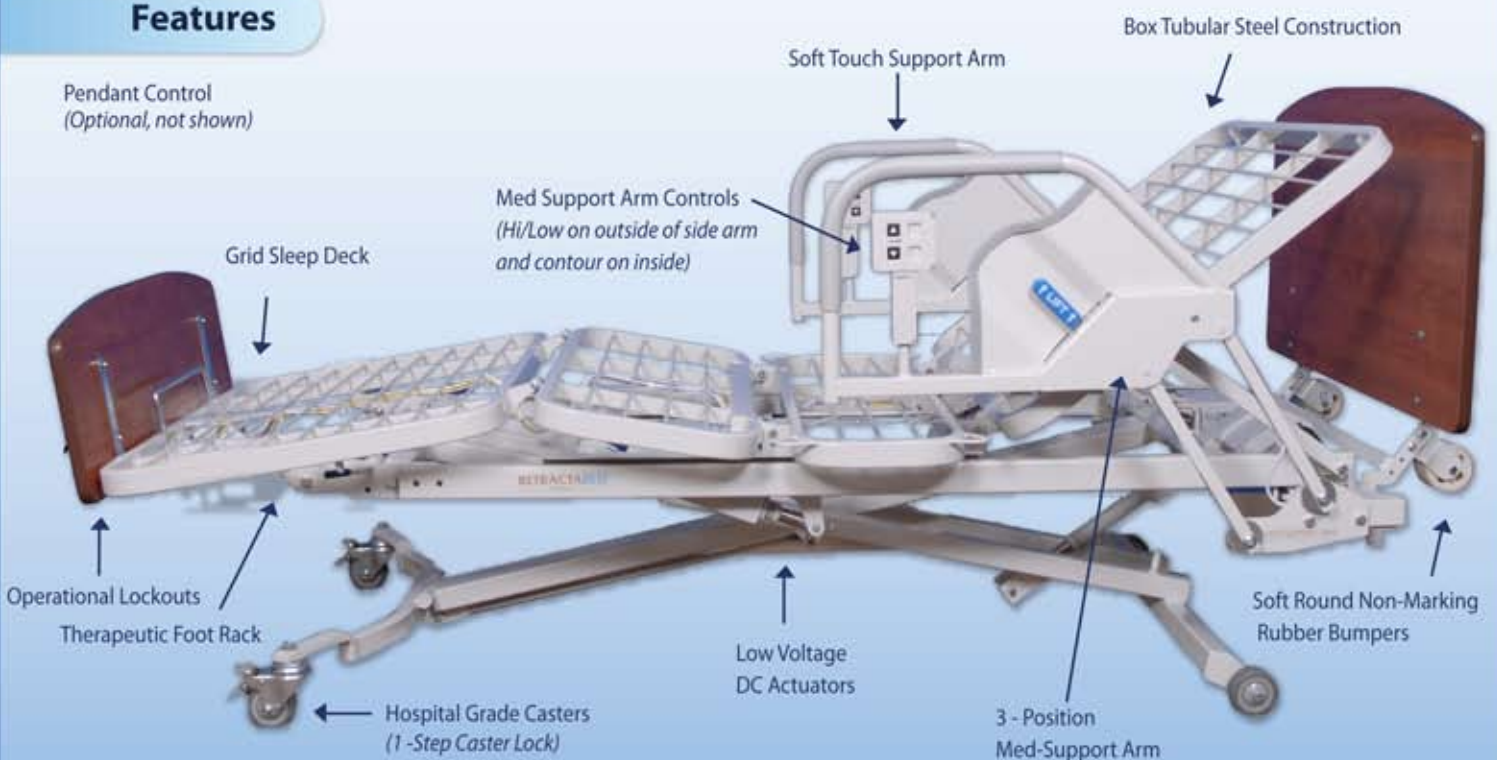

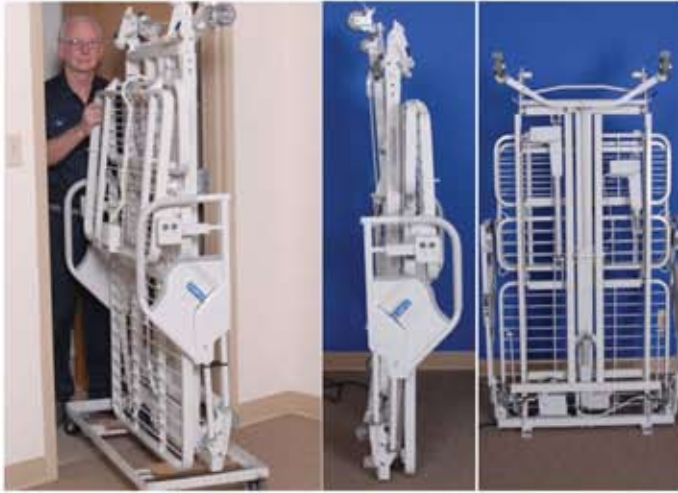

## Storage & Delivery Made Easy

**Setup:** When delivered by Med-Mizer, the RetractaBed® rolls off easily on four wheel carts and through any facility's doors or tight spaces. No tools are required for set-up and no cardboard boxes to discard.

**Storage:** Space is at a premium in any modern health care facility. The unique, space-saving Med-Mizer Stand-Up Storage feature permits a vertical footprint of approximately 4 square feet.

### Headboards to Match Any Decor



English Oak      Mahogany      Cherry      Light Oak

Med-Mizer offers standard head and foot boards in four different stains. However, custom head and foot boards are available.

### Specifications

	SS803A/803	SS763A/763
Deck Height Low	8.5"	8.5"
Deck Height High	24.5"	24.5"
Sleeping Surface	35" w x 80" L	35" w x 76" L
Length Retracted	71"	67"
Length Extended	85"	81"
Maximum Weight	450lbs	450lbs

## Quality Tested and Fully Assembled To Save You Money



Because Med-Mizer, Inc. is committed to providing a safe, easy to operate, and quality-built product, all of our beds are shipped fully assembled after having been quality inspected as an assembled unit. There is no need to use valuable time assembling the beds or disposing of shipping materials. This means that you can immediately place the beds in use after simply unbolting them from the delivery dolly/pallet and placing the head/footboards in place. Med-Mizer, Inc.'s assembly, testing and shipping process helps to save your facility valuable time and money and ensures that your residents have access to the beds when they need them.

### Find Out More Today

Call Toll Free: **1-877-867-7365**

Or Email: **sales@med-mizer.com**

### Warranty

15 year limited warranty on bed frame

3 years on drive system

1 year on wood products and casters

Lifetime warranty on all welds

**www.Med-Mizer.com**